
Rapid Detection of Carbapenemase-producing *Enterobacteriaceae*

Patrice Nordmann, Laurent Poirel,
and Laurent Dortet

To rapidly identify carbapenemase producers in *Enterobacteriaceae*, we developed the Carba NP test. The test uses isolated bacterial colonies and is based on in vitro hydrolysis of a carbapenem, imipenem. It was 100% sensitive and specific compared with molecular-based techniques. This rapid (<2 hours), inexpensive technique may be implemented in any laboratory.

Multidrug resistance is emerging worldwide at an alarming rate among a variety of bacterial species, causing both community-acquired and nosocomial infections (1). Carbapenems, the last line of therapy, are now frequently needed to treat nosocomial infections, and increasing resistance to this class of β -lactams leaves the health care system with almost no effective drugs (1). However, reports of carbapenem-resistant *Enterobacteriaceae* have increased (2,3). Resistance may be related to association of a decrease in bacterial outer-membrane permeability, with overexpression of β -lactamases with no carbapenemase activity or to expression of carbapenemases (2,4,5). Spread of carbapenemase producers is a relevant clinical issue because carbapenemases confer resistance to most β -lactams (2). Various carbapenemases have been reported in *Enterobacteriaceae*, such as the following types: *Klebsiella pneumoniae* carbapenemase (KPC; Ambler class A); Verona integron-encoded metallo- β -lactamase (VIM), imipenemase (IMP), New Delhi metallo- β -lactamase (NDM) (all Ambler class B); and oxacillinase-48 (OXA-48; Ambler class D) (2,4–6). In addition, carbapenemase producers are usually associated with many other non- β -lactam resistance determinants, which give rise to multidrug- and pandrug-resistant isolates (2,3,7).

Potential carbapenemase producers are currently screened first by susceptibility testing, using breakpoint values for carbapenems (2,8). However, this technique is time-consuming, and many carbapenemase producers do

not confer obvious resistance levels to carbapenems. There is a need for laboratories to search for carbapenemase producers (9). Phenotype-based techniques for identifying in vitro production of carbapenemase, such as the modified Hodge test, are not highly sensitive and specific (2,8,10). Detection of metallo- β -lactamase producers (IMP, VIM, NDM) and of KPC producers may be based on the inhibitory properties of several molecules but requires additional expertise and time (usually an extra 24–48 hours) (2,8,11,12). Furthermore, no inhibitors are available for detecting OXA-48-type producers that are spreading rapidly, at least in northern Africa, the Middle East, and Europe (2). Molecular detection of carbapenemase genes remains costly and requires substantial expertise. Both the phenotype-based techniques and molecular tests are time-consuming (at least 12–24 hours) and are poorly adapted to the clinical need for isolating patients rapidly to prevent nosocomial outbreaks.

We developed a novel test, described here, based on a technique designed to identify the hydrolysis of the β -lactam ring of a carbapenem. This test is rapid, sensitive and specific, and adaptable to any laboratory in the world.

The Study

We included in the study 162 carbapenemase-producing strains of various enterobacterial species isolated from clinical samples (e.g., blood cultures, urine, sputum) and of global origin (Table 1). This collection of strains also included 46 strains that were fully susceptible to carbapenems or showed a decreased susceptibility to carbapenems as a consequence of non-carbapenemase-based mechanisms (Table 2). Antibigrams were carried out for all strains on Mueller-Hinton agar (Becton Dickinson, Le Point de Chaix, France) according to guidelines of the Clinical and Laboratory Standards Institute (13). The Carba NP (Carbapenemase Nordmann-Poirel) test was performed as follows. One calibrated loop (10 μ L) of the tested strain directly recovered from the antibiogram was resuspended in a Tris-HCl 20 mmol/L lysis buffer (B-PERII, Bacterial Protein Extraction Reagent; Thermo Scientific Pierce, Rockford, IL, USA), vortexed for 1 minute and further incubated at room temperature for 30 minutes. This bacterial suspension was centrifuged at $10,000 \times g$ at room temperature for 5 minutes. Thirty μ L of the supernatant, corresponding to the enzymatic bacterial suspension, was mixed in a 96-well tray with 100 μ L of a 1-mL solution made of 3 mg of imipenem monohydrate (Sigma, Saint-Quentin Fallavier, France), pH 7.8, phenol red solution, and 0.1 mmol/L $ZnSO_4$ (Merck Millipore, Guyancourt, France). The phenol red solution was prepared by mixing 2 mL of a phenol red (Merck Millipore) solution 0.5% (wt/vol) with 16.6 mL of distilled water. The pH value was then adjusted to 7.8 by adding drops of 1 N NaOH. A mixture of the phe-

Author affiliations: Faculté de Médecine Paris Sud, Le Kremlin-Bicêtre, France (P. Nordmann, L. Poirel, L. Dortet); and Institut National de la Santé et de la Recherche Médicale, Paris, France (P. Nordmann, L. Poirel, L. Dortet)

DOI: <http://dx.doi.org/10.3201/eid1809.120355>

Table 1. Carbapenemase-producing clinical enterobacterial isolates subjected to the Carba NP test*

Ambler class, carbapenemase type	Species	β-Lactamase	No.	MIC range, mg/L			Carba NP test result	
				IMP	ERT	MER		
Class A								
KPC-type	<i>Klebsiella pneumoniae</i>	KPC-2	27	0.5->32	4->32	1->32	+	
		KPC-3	3	0.5-8	4->32	1-8	+	
	<i>Klebsiella ozaenae</i>	KPC-3	1	>32	>32	2	+	
	<i>Escherichia coli</i>	KPC-2	5	0.5-4	0.5>32	0.5-2	+	
	<i>Enterobacter cloacae</i>	KPC-2	7	1-24	1.5-32	0.75-16	+	
	<i>Enterobacter aerogenes</i>	KPC-2	1	8	>32	8	+	
	<i>Citrobacter freundii</i>	KPC-2	2	8->32	1.5->32	1.5-3	+	
	<i>Serratia marcescens</i>	KPC-2	2	>32	>32	>32	+	
	<i>Salmonella</i> spp.	KPC-2	1	4	1	0.25	+	
	NMC-A	<i>E. cloacae</i>	NMC-A	1	16	>32	16	+
	SME-type	<i>S. marcescens</i>	SME-1	1	32	4	12	+
			SME-2	1	32	4	12	+
	GES-type	<i>E. cloacae</i>	GES-5	1	>32	>32	>32	+
IMI-type	<i>Enterobacter asburiae</i>	IMI-2	1	>32	>32	>32	+	
Class B								
NDM-type	<i>K. pneumoniae</i>	NDM-1	16	0.5->32	2->32	1->32	+	
		NDM-4	1	>32	>32	>32	+	
	<i>E. coli</i>	NDM-1	7	1-16	3->32	1-16	+	
	<i>E. cloacae</i>	NDM-1	1	2	16	2	+	
	<i>C. freundii</i>	NDM-1	1	>32	>32	>32	+	
	<i>Providencia stuartii</i>	NDM-1	1	12	0.38	1.5	+	
	<i>Proteus rettgeri</i>	NDM-1	1	3	0.5	1.5	+	
	VIM-type	<i>K. pneumoniae</i>	VIM-1	15	0.5->32	0.5->32	0.38- >32	+
			VIM-19	1	8	16	4	+
		<i>E. coli</i>	VIM-1	2	1.5-3	0.38-1.5	0.5-1	+
VIM-2			2	2-4	0.5-1.5	0.38- 0.5	+	
<i>E. cloacae</i>		VIM-19	1	8	16	4	+	
		VIM-1	4	1->32	0.38 to >32	0.5->32	+	
IMP-type	<i>S. marcescens</i>	VIM-2	1	>32	>32	>32	+	
	<i>K. pneumoniae</i>	IMP-1	5	0.5-8	2-4	1-8	+	
		IMP-8	2	0.5-1	0.5-1	0.5	+	
	<i>E. coli</i>	IMP-1	2	0.5	3-4	0.5-1	+	
		IMP-8	1	6	8	3	+	
	<i>E. cloacae</i>	IMP-1	12	8->32	>32	2->32	+	
		IMP-8	2	0.75- 1.5	0.5-1	0.5-1	+	
	<i>S. marcescens</i>	IMP-1	2	8->32	>32	2->32	+	
		IMP-11	1	8	>32	2	+	
	Class D							
OXA-48 type	<i>K. pneumoniae</i>	OXA-48	15	0.38- >32	0.38->32	0.38- >32	+	
		OXA-181	2	0.5-1	2-4	0.5-1	+	
	<i>E. coli</i>	OXA-48	6	0.38-3	0.5-16	0.12-1	+	
	<i>E. cloacae</i>	OXA-48	3	0.5-1	0.5-16	0.5-1.5	+	
	<i>P. rettgeri</i>	OXA-181	1	8	1	2	+	

*IMP, imipenem; ERT, ertapenem; MER, meropenem; KPC, *K. pneumoniae* carbapenemase; NMC-A, non-metallo-enzyme carbapenemase; SME, *S. marcescens* enzyme; GES, Guiana extended-spectrum β-lactamase; IMI, imipenem-hydrolysing β-lactamase; NDM, New Delhi metallo-β-lactamase; VIM, Verona integron-encoded metallo-β-lactamase; IMP, imipenemase; OXA-48, oxacillinase-48.

nol red solution and the enzymatic suspension being tested was incubated at 37°C for a maximum of 2 hours. Test results were interpreted by technicians who were blinded to the identity of the patients who gave the samples.

All strains had previously been characterized for their β-lactamase content at the molecular level. MICs of carbapenems were determined by using the Etest (AB bioMérieux, Solna, Sweden), and results were recorded according to US guidelines (Clinical and Laboratory Standards Institute), as updated in 2012 (13). The breakpoints used were

those for imipenem and meropenem: susceptibility <1 μg/mL, resistance >4 μg/mL, and for ertapenem, susceptibility <0.25 μg/mL, resistance > μg/mL.

When the Carba NP test was used, the color of the wells turned from red to orange or yellow (Figure 1) for all tested strains that were producing carbapenemases (Table 1), whereas wells corresponding to bacterial extracts of isolates that did not produce carbapenemase remained red, whatever their level of carbapenem susceptibility (Table 2). The color changed from red to yellow as early as 5–10

Table 2. Non-carbapenemase-producing clinical enterobacterial isolates subjected to the Carba NP test*

β-Lactamase type, species	β-Lactamase	No.	MIC, mg/L			Carba NP test result
			IMP	ERT	MER	
ESBLs						
<i>Klebsiella pneumoniae</i>	CTX-M-3	1	0.12	0.12	0.12	–
	CTX-M-14	1	0.12	0.12	0.12	–
	CTX-M-15	3	0.12	0.12	0.12	–
<i>Escherichia coli</i>	CTX-M-1	1	0.12	0.12	0.12	–
	CTX-M-3	1	0.12	0.12	0.12	–
	CTX-M-14	2	0.12	0.12	0.12	–
	CTX-M-15	2	0.12	0.12	0.12	–
	VEB-1	1	0.12–0.25	0.12	0.12	–
<i>Enterobacter cloacae</i>	CTX-M-15	3	0.12	0.12	0.12	–
	VEB-1	1	0.12	0.12	0.12	–
Plasmid-mediated AmpC or chromosomal AmpC + decreased membrane permeability						
<i>K. pneumoniae</i>	DHA-1	1	>32	>32	>32	–
	DHA-2	1	0.12	0.5	0.12	–
<i>E. coli</i>	Extended-spectrum cephalosporinase	1	0.12	0.12	0.12	–
	CMY-2	1	0.12	0.12	0.12	–
	CMY-10	1	0.12	0.38	0.12	–
	DHA-1	1	0.12	0.12	0.12	–
	ACC-1	1	0.12	0.12	0.12	–
	Overexpressed cephalosporinase	1	16	>32	2	–
<i>Proteus mirabilis</i>	ACC-1	1	0.25	0.12	0.12	–
<i>E. cloacae</i>	Overexpressed cephalosporinase	7	0.12–16	1–>32	0.12–>32	–
<i>Enterobacter aerogenes</i>	Overexpressed cephalosporinase	1	1	4	0.75	–
<i>Morganella morganii</i>	Overexpressed cephalosporinase	2	1.5–2	0.12	0.5	–
ESBL + decreased membrane permeability						
<i>K. pneumoniae</i>	CTX-M-15	8	0.25–8	1–>32	1–>32	–
	SHV-28	1	1	4	1	–
	SHV-2a	1	0.25	2	0.38	–
<i>Enterobacter sakazakii</i>	CTX-M-15	1	0.25	1.5	0.25	–
<i>Citrobacter freundii</i>	TEM-3	1	1	8	1	–

*IMP, imipenem; ERT, ertapenem; MER, meropenem; ESBLs, extended-spectrum β-lactamases.

minutes after incubation for KPC producers began. In most cases, incubation for 30 minutes was sufficient for obtaining a frank color change for carbapenemase producers. The test's specificity and sensitivity were 100% when results

were compared with those from molecular-based methods, the reference standard for identifying carbapenemase genes. All tests were performed in triplicate, giving identical and reproducible results.

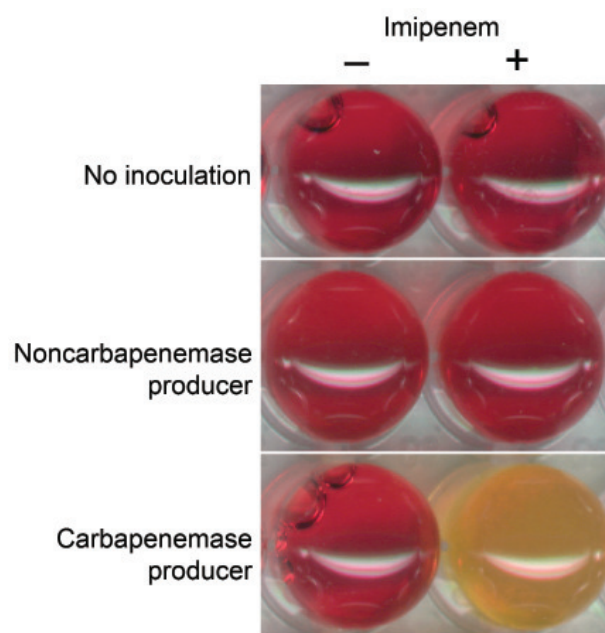


Figure 1. Representative results of the Carba NP test. The Carba NP test was performed with a noncarbapenemase producer (*Escherichia coli* producing the extended-spectrum β-lactamase CTX-M-15, upper panel) and with a carbapenemase producer (*Klebsiella pneumoniae*–producing New Delhi metallo-β-lactamase-1, lower panel) in a reaction medium without (left panel) and with (right panel) imipenem. Uninoculated wells are shown as controls. Photographs were taken after a 1.5-hour incubation. A color version of this figure is available online (wwwnc.cdc.gov/EID/article/18/9/12-0355-F1.htm).

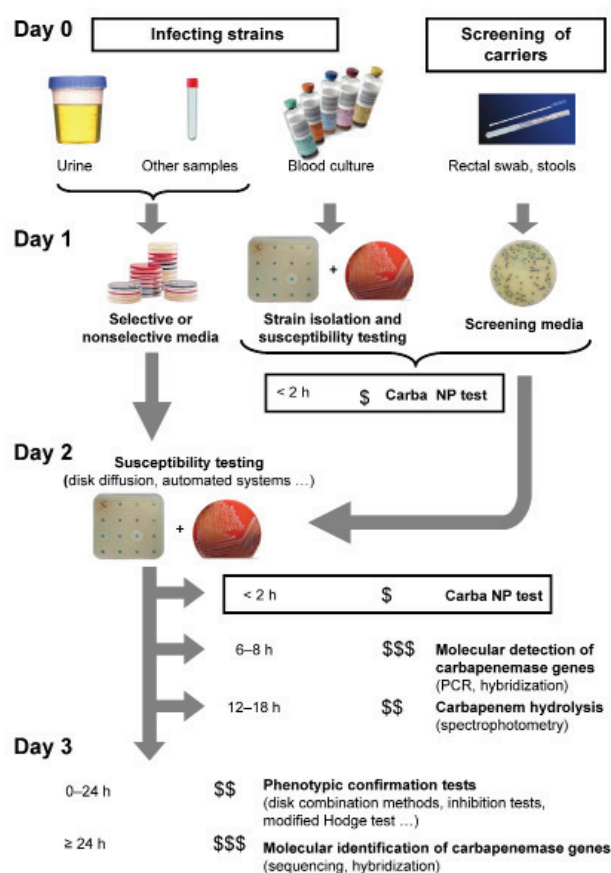


Figure 2. Strategy for identification of carbapenemase-producing *Enterobacteriaceae*. The time needed to perform the test is indicated before each test. The number of flasks indicates the degree of specialization needed to perform the test; the number of \$ indicates the relative cost of each test. A color version of this figure is available online (wwwnc.cdc.gov/EID/article/18/9/12-0355-F2.htm).

The Carba NP test perfectly differentiates carbapenemase producers (Table 1) from strains that are carbapenem resistant due to non-carbapenem-mediated mechanisms, such as combined mechanisms of resistance (outer-membrane permeability defect associated with overproduction of cephalosporinase and/or extended-spectrum β -lactamases) or from strains that are carbapenem susceptible but express a broad-spectrum β -lactamase without carbapenemase activity (extended-spectrum β -lactamase, plasmid and chromosome-encoded cephalosporinases) (Table 2). Interpretable positive results were obtained in <2 hours, making it possible to implement rapid containment measures to limit the spread of carbapenemase producers.

Conclusions

The Carba NP test has multiple benefits. It is inexpensive, rapid, reproducible, and highly sensitive and specific.

It eliminates the need for using other techniques to identify carbapenemase producers that are time-consuming and less sensitive or specific. Using this accurate test would improve detection of patients infected or colonized with carbapenemase producers. The test has been routinely implemented in our microbiology department at Hôpital de Bicêtre and is giving excellent results (data not shown). In addition, use of the Carba NP test has greatly decreased the laboratory technicians' workload and simplified the clinical management of potential carbapenemase producers.

This test could be used, for example, for directly testing 1) bacteria obtained from antibiograms of blood culture or 2) bacterial colonies grown on culture media before antimicrobial drug susceptibility testing (Figure 2). Further studies will evaluate its clinical value for antimicrobial drug stewardship on bacteria isolated directly from clinical samples (Figure 2). When the Carba NP test is used for that purpose, we expect that the time to detect carbapenemase producers will decrease by at least 24 hours (Figure 2).

The test could also be used to quickly identify carbapenem-resistant isolates from fecal specimens screened for multidrug-resistant bacteria (Figure 2). This capability would be valuable in preventing outbreaks. To determine positive and negative predictive values of the test, additional evaluations will be required with strains isolated from clinical samples screened on different types of selective media. The use of the Carba NP test may also support novel antimicrobial drug development by facilitating patient enrollment in pivotal clinical trials. Its use as a home-made test may contribute to the global surveillance network.

The Carba NP test can efficiently indicate the strains to be further tested by PCR or submitted to sequencing for a detailed identification of the carbapenemase genes. Last, the test could be used in low-income countries that are large reservoirs for carbapenemase producers (2). It offers a practical solution for detecting a main component of multidrug resistance in *Enterobacteriaceae*. Use of the Carba NP test will contribute to a better stewardship of carbapenems by changing the paradigm of controlling carbapenemase producers worldwide

This work was funded by a grant from the Institut National de la Santé et de la Recherche Médicale (INSERM) (U914). An international patent form has been filed on behalf of INSERM Transfert (Paris, France) corresponding to this test.

Dr Nordmann is professor of medical microbiology, South-Paris University, Paris, and director of the INSERM U914 Emerging Resistance to Antibiotics program. His main field of research interest includes the genetics, biochemistry, and molecular epidemiology of resistance in gram-negative bacteria.

References

1. Infectious Diseases Society of America; Spellberg B, Blaser M, Guidos RJ, Boucher HW, Bradley JS, et al. Combating antimicrobial resistance: policy recommendations to save lives. *Clin Infect Dis*. 2011;52(Suppl 5):S397–428. <http://dx.doi.org/10.1093/cid/cir153>
2. Nordmann P, Naas T, Poirel L. Global spread of carbapenemase-producing *Enterobacteriaceae*. *Emerg Infect Dis*. 2011;17:1791–8.
3. Schwaber MJ, Carmeli Y. Carbapenem-resistant *Enterobacteriaceae*: a potential threat. *JAMA*. 2008;300:2911–3. <http://dx.doi.org/10.1001/jama.2008.896>
4. Martínez-Martínez L. Extended-spectrum β -lactamases and the permeability barrier. *Clin Microbiol Infect*. 2008;14(Suppl 1):82–9. <http://dx.doi.org/10.1111/j.1469-0691.2007.01860.x>
5. Thomson KS. Extended-spectrum β -lactamase, AmpC, and carbapenemase issues. *J Clin Microbiol*. 2010;48:1019–25. <http://dx.doi.org/10.1128/JCM.00219-10>
6. Moellering RC Jr. NDM-1—a cause for worldwide concern. *N Engl J Med*. 2010;363:2377–9. <http://dx.doi.org/10.1056/NEJMp1011715>
7. Walsh TR, Toleman MA. The emergence of pan-resistant gram-negative pathogens merits a rapid global political response. *J Antimicrob Chemother*. 2012;67:1–3. <http://dx.doi.org/10.1093/jac/dkr378>
8. Nordmann P, Gniadkowski M, Giske CG, Poirel L, Woodford N, Miriagou V. European Network on Carbapenemases. Identification and screening of carbapenemase-producing *Enterobacteriaceae*. *Clin Microbiol Infect*. 2012;18:432–8. <http://dx.doi.org/10.1111/j.1469-0691.2012.03815.x>
9. Livermore DM, Andrews JM, Hawkey PM, Ho PL, Keness Y, Doi Y, et al. Are susceptibility tests enough, or should laboratories still seek ESBLs and carbapenemases directly? *J Antimicrob Chemother*. 2012;67:1569–77. <http://dx.doi.org/10.1093/jac/dks088>
10. Pasteran F, Mendez T, Rapoport M, Guerriero L, Corso A. Controlling false-positive results obtained with the Hodge and Masuda assays for detection of class A carbapenemase in species of *Enterobacteriaceae* by incorporating boronic acid. *J Clin Microbiol*. 2010;48:1323–32. <http://dx.doi.org/10.1128/JCM.01771-09>
11. Centers of Disease Control and Prevention. Detection of *Enterobacteriaceae* isolates carrying metallo- β -lactamase—United States, 2010. *MMWR Morb Mortal Wkly Rep*. 2010;59:750.
12. Nordmann P, Poirel L, Carrè A, Toleman MA, Walsh TR. How to detect NDM-1 producers? *J Clin Microbiol*. 2011;49:718–21. <http://dx.doi.org/10.1128/JCM.01773-10>
13. Clinical and Laboratory Standards Institute. Performance standards for antimicrobial susceptibility testing. M100–S22. Wayne (PA): the Institute; 2012.

Address for correspondence: Patrice Nordmann, Service de Bactériologie-Virologie, Hôpital de Bicêtre 78 Rue du Général Leclerc, 94275 K-Bicêtre, France; email: nordmann.patrice@bct.aphp.fr

Health-e-Cards: Welcome - Windows Internet Explorer

<http://www2.cdc.gov/ecards/>

Send your colleagues, family, and friends eCards so they can find out about the latest emerging infectious diseases

CDC has more than 100 free Health-e-Cards (or "electronic greeting cards") to send to friends, family, and co-workers! Simply select a CDC Health-e-Card in an email address and name, and send a colorful greeting to promote healthy living, promote safe activities, or celebrate a special event!

Popular eCards:

- Diseases
- Environment
- Holidays
- Travel

New eCards:

- Keep Teens Healthy
- Blooming Baby Healthy Living
- Voices are Blue
- Happy Valentine's Day from Healthcare.gov
- EID Journal Diseases & Conditions
- Discover the Icy Realm of the Rime
- Fruit and Vegetable Budget Tips
- March is National Nutrition Month. Send an eCard to learn how you can save money and eat healthy.

Text size: [A] [B] [C] [D] [E] [F] [G] [H] [I] [J] [K] [L] [M] [N] [O] [P] [Q] [R] [S] [T] [U] [V] [W] [X] [Y] [Z]

Email page
Print page
Bookmark and Share

View page in:
Español (Spanish)

Contact Us:
Centers for Disease Control and Prevention
1600 Clifton Rd
Atlanta, GA 30333
800-CDC-INFO (800-232-4636)
TTY: (888) 232-6348
24 Hours/Every Day
cdcinfo@cdc.gov

How are we doing?
Tell us what you think about Health-e-Cards >>

Trusted sites 75%